

This Page Is Inserted by IFW Operations
and is not a part of the Official Record

BEST AVAILABLE IMAGES

Defective images within this document are accurate representations of the original documents submitted by the applicant.

Defects in the images may include (but are not limited to):

- BLACK BORDERS
- TEXT CUT OFF AT TOP, BOTTOM OR SIDES
- FADED TEXT
- ILLEGIBLE TEXT
- SKEWED/SLANTED IMAGES
- COLORED PHOTOS
- BLACK OR VERY BLACK AND WHITE DARK PHOTOS
- GRAY SCALE DOCUMENTS

IMAGES ARE BEST AVAILABLE COPY.

**As rescanning documents *will not* correct images,
please do not report the images to the
Image Problem Mailbox.**

Bl66



PCT

WORLD INTELLECTUAL PROPERTY ORGANIZATION
International Bureau

INTERNATIONAL APPLICATION PUBLISHED UNDER THE PATENT COOPERATION TREATY. (PCT)

(51) International Patent Classification ⁶ : C12Q 1/68		A1	(11) International Publication Number: WO 99/20798 (43) International Publication Date: 29 April 1999 (29.04.99)
(21) International Application Number: PCT/US98/22406 (22) International Filing Date: 22 October 1998 (22.10.98) (30) Priority Data: 60/063,219 23 October 1997 (23.10.97) US (71) Applicant: EXACT LABORATORIES, INC. [US/US]; 63 Great Road, Maynard, MA 01754 (US). (72) Inventor: SHUBER, Anthony, P.; 11 Grant Street, Milford, MA 01757 (US). (74) Agent: LINKKILA, Timothy, P.; Testa, Hurwitz & Thibault, LLP, High Street Tower, 125 High Street, Boston, MA 02110 (US).			(81) Designated States: AU, CA, JP, European patent (AT, BE, CH, CY, DE, DK, ES, FI, FR, GB, GR, IE, IT, LU, MC, NL, PT, SE). Published <i>With international search report. Before the expiration of the time limit for amending the claims and to be republished in the event of the receipt of amendments.</i>
(54) Title: METHODS FOR DETECTING CONTAMINATION IN MOLECULAR DIAGNOSTICS USING PCR			
<p>The diagram illustrates the PCR process. At the top, a template DNA strand is shown with a 5' end on the left and a 3' end on the right. A 'Forward, ss' primer is shown binding to the 5' end of the template, and a 'Rev., ss' primer is shown binding to the 3' end. The template is labeled 'CDS'. Below the template, a vertical arrow labeled 'PCR' points to the resulting PCR products. The PCR products are shown as three horizontal lines, each with a 'CDS' label at both ends, indicating the amplified DNA fragments.</p>			
(57) Abstract			
<p>The invention provides methods for detecting contamination in a PCR reaction. Methods of the invention are especially useful for detection of contamination in heterogeneous samples containing a rare nucleic acid to be detected.</p>			

FOR THE PURPOSES OF INFORMATION ONLY

Codes used to identify States party to the PCT on the front pages of pamphlets publishing international applications under the PCT.

AL	Albania	ES	Spain	LS	Lesotho	SI	Slovenia
AM	Armenia	FI	Finland	LT	Lithuania	SK	Slovakia
AT	Austria	FR	France	LU	Luxembourg	SN	Senegal
AU	Australia	GA	Gabon	LV	Latvia	SZ	Swaziland
AZ	Azerbaijan	GB	United Kingdom	MC	Monaco	TD	Chad
BA	Bosnia and Herzegovina	GE	Georgia	MD	Republic of Moldova	TG	Togo
BB	Barbados	GH	Ghana	MG	Madagascar	TJ	Tajikistan
BE	Belgium	GN	Guinea	MK	The former Yugoslav Republic of Macedonia	TM	Turkmenistan
BF	Burkina Faso	GR	Greece	ML	Mali	TR	Turkey
BG	Bulgaria	HU	Hungary	MN	Mongolia	TT	Trinidad and Tobago
BJ	Benin	IE	Ireland	MR	Mauritania	UA	Ukraine
BR	Brazil	IL	Israel	MW	Malawi	UG	Uganda
BY	Belarus	IS	Iceland	MX	Mexico	US	United States of America
CA	Canada	IT	Italy	NE	Niger	UZ	Uzbekistan
CF	Central African Republic	JP	Japan	NL	Netherlands	VN	Viet Nam
CG	Congo	KE	Kenya	NO	Norway	YU	Yugoslavia
CH	Switzerland	KG	Kyrgyzstan	NZ	New Zealand	ZW	Zimbabwe
CI	Côte d'Ivoire	KP	Democratic People's Republic of Korea	PL	Poland		
CM	Cameroon	KR	Republic of Korea	PT	Portugal		
CN	China	KZ	Kazakhstan	RO	Romania		
CU	Cuba	LC	Saint Lucia	RU	Russian Federation		
CZ	Czech Republic	LI	Liechtenstein	SD	Sudan		
DE	Germany	LK	Sri Lanka	SE	Sweden		
DK	Denmark	LR	Liberia	SG	Singapore		
EE	Estonia						

METHODS FOR DETECTING CONTAMINATION IN MOLECULAR DIAGNOSTICS USING PCR

Background of the Invention

The polymerase chain reaction (PCR) is a widely-used clinical laboratory procedure for sequence-specific target amplification. However, contamination is an ongoing problem. For many PCR applications, it is essential that the only DNA that enters the reaction is the template to
5 be amplified.

Increases in the sensitivity and specificity of PCR have enabled analysis of heterogeneous DNA (e.g., from tumor biopsies, stool). The DNA to be amplified is typically a rare event in the context of a heterogeneous sample. However, as the degree of sample heterogeneity increases, the tolerable threshold of background (signal generated from a negative control sample) becomes
10 increasingly lower. This is necessary to retain a sufficient signal to noise ratio between positive clinical samples and negative control samples within an assay, and therefore to retain high confidence in the assay results. The end result of applying PCR to more heterogeneous DNA environments is a reduced tolerance for pre-PCR contamination from previous amplified material. Currently, there are three methods applied to prevent PCR contamination: (1) physical separation
15 of the sample, pre-PCR setup, and post-PCR manipulations; (2) use of Uracil DNA-glycosylase and dUTP instead of dTTP, and (3) the use of UV irradiation.

Thousands of samples may be analyzed in a single clinical assay with multiple PCR negative controls added. In this context, an investigator relies on the presence or absence of amplified product within a limited number of negative control samples to confirm the origin of
20 amplification products observed in experimental samples. If only one PCR negative control sample is positive, the entire assay is invalid, and must be repeated. In an assay containing 1000 samples, each sample must be run with another set of negative controls when contamination is observed.

However, the mere lack of amplification product within the PCR negative control is not
25 determinative of a positive PCR result in a sample in which contamination is rare. This kind of sporadic contamination is especially problematic in an extremely large throughput assay in which 5 to 10 negative controls are run for approximately every 1000 samples. Statistically, the likelihood of sporadic contamination in, for example, 1000 samples will not be detected in only 5

- 2 -

negative controls. Sporadic contamination is also a significant problem when PCR based analyses are performed on heterogeneous (rare event analysis) samples in which a positive result is generated from, for example, 1-5% of the total amplification product present within the sample. Generally, within a PCR based inherited disease diagnostic assay, given the 50% heterogeneity that exists in any genomic DNA sample, a 1-5% increase in signal in a true negative sample would appear as a slight increase in background, but would not indicate a false positive result. However, within an assay involving samples with heterogeneous populations of DNA, a 1-5% positive signal generated by a true negative sample would result in a false positive.

In addition, even within an inherited disease diagnostic assay, if there were 1000 samples analyzed and 5-10 negative control PCR reactions were run in parallel, and one or two of the negative control samples were positive, results from any of the samples themselves would be compromised. If the contamination of the PCR negative control samples is truly sporadic, then repeat analysis of all 1000 samples is probably not necessary and extremely costly. The lack of amplification product within the PCR negative control samples is not determinative that a positive PCR result within an experimental sample set is not from rare (sporadic) contamination that has occurred in only a few samples within the assay (and not due to the negative controls run in parallel).

In many assays, "normal" PCR contaminants (e.g., resulting from purification problems) are an even greater hindrance and leads to decreased sensitivity of the assay. These "normal" PCR contaminants can lead to false negative results that undermine the accuracy of (and confidence in) the particular assay.

Therefore, methods are needed for performing clinical analyses on samples of DNA heterogeneity (e.g. sporadic cancer detection) such that sporadic contamination from previous amplification product or "normal" PCR contaminants do not result in false positive or false negative results.

Summary of the Invention

The invention provides methods for determining whether contamination from previous amplification product exists in products of a polymerase chain reaction (PCR). Specifically, the invention relates to methods for detecting the presence of PCR products (amplicons) that would not be present *but for* contamination from previous amplification product in the PCR sample. Methods of the invention are useful for detection of contamination in any PCR. Furthermore, the

- 3 -

methods of the invention are useful to avoid false negative and false positive results and the decreased assay sensitivity associated with PCR contamination. However, such methods are especially useful in heterogeneous samples, particularly samples in which the detection of a rare event (i.e. a small subpopulation of a nucleic acid in a heterogeneous sample) is the ultimate
5 object of the PCR.

In a preferred embodiment, methods of the invention comprise the utilization of optimal primer construction for PCR. Accordingly, in a highly-preferred embodiment, methods of the invention comprise conducting a first amplification using one or more (preferably two) chimeric primers. A chimeric primer, for purposes of the invention, is one comprising a primer having
10 substantial sequence specificity with the template to be amplified (a template-specific sequence) and a 5' end that is referred to herein as a "contamination detection sequence" (CDS). Methods further comprise conducting a second, parallel, amplification reaction using at least one (preferably two) contamination detection sequence (without the attached template-specific sequence, or with only a minimal number of template-specific bases, as described below) as a
15 primer. Finally, detection of an amplicon in the second amplification reaction means that the sample is contaminated with previous amplification product, because no such amplicon would be generated absent contamination.

A contamination detection sequence may be any sequence (regardless of length) that does not have substantial sequence specificity (i.e., does not hybridize under stringent conditions) with
20 the template. See Figure 2A. Accordingly, methods of the invention comprise conducting two amplification reactions on each sample suspected to contain a template sequence, the amplification of which is desired. The first reaction utilizes the chimeric primers described above, and yields the desired amplicon (which may then be sequenced, probed, etc.). The second amplification reaction utilizes only the contamination detection sequence as primers (which are
25 non-specific relative to the template). Any amplicon produced in the second reaction is evidence of contamination with previous amplification product.

In a preferred embodiment, the amplification reaction is selected from PCR, reverse transcriptase PCR, and Q-PCR. Also in a preferred embodiment, the sample containing nucleic acid to be amplified is a stool sample. A stool sample contains a highly-heterogeneous population
30 of nucleic acids. Human nucleic acids represent a small portion of the nucleic acid present in stool. More specifically, a stool sample may contain molecular indicia of cancer, specifically

- 4 -

colorectal cancer, that occurs as a small subpopulation (typically on the order of about 1% at early stages of cancer or precancer) of the total nucleic acid in the stool. Sensitive assays (which may or may not involve amplification) have been developed to detect such small subpopulations. See, e.g., U.S. Patent No. 5,670,325, incorporated by reference herein. Amplification of a nucleic acid containing a mutation indicative of cancer or precancer may be confounded by PCR contaminants in the sample, especially if the detection limits of the assay are near or above the percent contaminants in the sample. The present invention detects PCR contaminants, thus allowing a given PCR reaction to be excluded from analysis on a sample-by-sample basis. Thus, if four separate samples are taken for amplification, each sample is divided into two subsample aliquots, one of which is amplified using chimeric primers, and in the other, the contamination detection sequence primers are used to check for contamination in the aliquot. Therefore, each aliquot of sample for which amplification is sought has its own quality control assay.

These and other advantages and aspects of the invention will be understood upon consideration of the following detailed description thereof.

15 **Brief Description of the Drawings**

Figure 1 shows a schematic diagram of the chimeric primers used in the present invention.

Figure 2A shows a schematic representation of PCR amplification using chimeric primers.

Figure 2B shows a schematic representation of PCR amplification using CDS-specific primers on uncontaminated sample (resulting in no amplified product).

20 Figure 2C shows a schematic representation of PCR amplification using CDS-specific primers on contaminated sample (resulting in amplified product).

Figure 3 shows a schematic diagram of an assay performed using the methods of the invention.

Detailed Description of the Invention

25 Methods of the invention comprise optimal PCR primer design. Normally, target-specific PCR primers are complementary to sequences present within the target. The target sequence is part of, and endogenous to, the target DNA analyte (the analyte can be any target DNA of interest: human DNA, viral DNA etc.) and is therefore one that is expected to be present in all of the experimental samples (in the case of inherited disease diagnostics), or at least in all of the positive samples (in the case of infectious disease diagnostics).

30

- 5 -

Methods of the invention comprise PCR primers that have a non-homologous or non-complementary "contamination detection sequence" ("CDS") attached to the 5' end of target-specific PCR primers. (See Fig. 1). The CDS region is neither homologous to, nor complementary to, any endogenase (template) sequence. Therefore, following any PCR involving chimeric primers, the CDS becomes incorporated into the PCR products (amplicons) generated from the PCR. (See Fig. 2A). Therefore, only PCR products from previous reactions have the CDS region contained within them.

In the present invention, a sample to be assayed for a particular analyte (which may be one of hundreds or thousands in a single clinical assay) is analyzed by two distinct, parallel amplification reactions. In a first reaction, PCR is performed on the sample using chimeric sequences that contain a template-specific sequence (a sequence substantially complementary to a specific DNA analyte) and a 5' CDS sequence (the CDS sequence is contiguous to the 5' end of the template specific sequence). In a second reaction, PCR is performed on the sample using primers that are specific for previously amplified amplicons containing the CDS sequence. The CDS primers (1) may be sequences that are specific for the CDS sequence alone (i.e., with no cross-reactivity to the target analyte sequence), or (2) the primers can comprise the CDS sequence with additional bases attached to the 3' end. From one to nine bases may be added at the 3' end of the CDS sequence and may serve to provide additional specificity. The CDS primers will not effectively prime the target analyte sequence.

The presence of amplified product (or amplicon) as a result of the PCR using the chimeric primer indicates a positive result for the presence of the particular analyte, but may also reflect contamination from previous PCR product.

The PCR with CDS-specific primers acts as negative control. Because the only samples that contain the CDS sequence will be those generated by previous PCR events within the lab, the presence of amplified product after PCR with the CDS-specific primers indicates that that particular sample is contaminated, and the results should be discarded. The lack of amplified product reflects the absence of PCR-based contaminants in a particular sample as portion of a sample.

Accordingly, the present invention eliminates the degree of repeat sample analyses performed within, for example, high throughput assays by specifically identifying only samples that have contamination.

- 6 -

The invention also provides additional protection needed within assays of heterogeneous samples, where sporadic contamination is more likely to be the source of contamination. It enables identification of the specific samples within the essays that are truly contaminated. The invention also provides a sample specific internal control for determining PCR product contamination.

The present invention is suitable for use with a variety of experimental samples that may contain a particular DNA analyte. Biological samples may be used in the present invention, including blood and stool samples.

The methods of the present invention are especially suitable for applications such as inherited disease diagnostics and related kits; infectious disease diagnostics and related kits; clinical assays involving sporadic cancer detection (e.g. testing DNA from stool for colorectal cancer) and related kits; and other "rare event" clinical assay and related kits.

Example 1

Stool sample is collected and prepared as described in U.S. Patent No. 5,741,650 and in copending applications Serial No. 08/876,638, filed June 16, 1997, both incorporated by reference herein. Specifically, stool is collected and prepared so that a sample contains at least a cross-sectional portion of a stool voided by a patient. Alternatively, whole stool may be used. The sample is homogenized in a physiologically compatible buffer (e.g., having a final concentration: 500 mM Tris, 16 mM EDTA and 10 mM NaCl, pH 9.0), using an Exactor II shaker for 15 minutes. A 20% SDS solution is added to a final concentration of 0.5%. Proteinase K is also added to a final concentration of 500 µg/ml and incubated at 37°C.

For exemplification, sequence-specific primers suitable for PCR are chosen to correspond to a portion of the *kras* gene sequence. These are: Primer 1 (SEQ. ID. NO. 2): 5'-GATTCCTACA GGAAGCAAGTAGTAATTG-3', and Primer 2 (SEQ. ID. NO. 3): 5'-TAATGGTGAATATCTTCAAATGATTTAG-3'.

The contamination detection sequence (CDS) is 5'-GCGGTCCCAAAGGGTCAGT-3' (SEQ. ID. NO. 1). The chimeric primers contain the 20-nucleotide CDS sequence attached (i.e. contiguous) to the 5' end of the individual sequence-specific primers (primer 1 or primer 2). Oligonucleotides are HPLC purified and quantitated by spectrophotometry.

PCR amplifications are performed using from about 4 µl (1-2 µg) to about 10 µl (5-50 ng) of genomic DNA prepared from stool samples. PCR amplifications are done using a Perkin-

- 7 -

Elmer 9600 Thermal Cycler (Perkin-Elmer, Norwalk, CT) for 28 cycles with ramping (94°C/10-sec hold with 48-sec ramp, 60°C/10-sec hold with 36-sec ramp, 72°C/10-sec hold with 38 sec ramp). Reactions (50 µl) are carried out in 1 x PCR buffer (10 mM Tris-HCl at pH 8.3, 50 mM KCl, 1.5 mM MgCl₂), 200 µM dNTPs, 2.5 units, of *Taq* polymerase (Perkin-Elmer, Norwalk, CT).

For PCR product analyses, 8 µl of the amplification reactions is loaded directly onto a 2% ethidium bromide stained agarose gel and electrophoresed at 250 V for 90 min. The amplification products are visualized with a UV transilluminator (Fotodyne, New Berlin, WI) and photographed with an Alpha Innotech IS-500 Digital Imaging System version 1.97 (Sun Bioscience Inc., Branford, CT).

A first PCR is performed in a first aliquot of stool sample (containing *kras*) using chimeric primers, wherein the forward primer comprises primer 1 with the CDS contiguous with its 5' end; and the reverse primer comprises primer 2 with the CDS contiguous with its 5' end. The first PCR results in an amplicon comprising both chimeric primers and the intervening template sequence.

A second PCR is performed on a second aliquot of stool sample in which both the forward and reverse primers are the CDS. If contamination from previous PCR cycles is present in the sample, the second PCR will product an amplicon. If no contamination is present in the sample, no amplicon is observed in the second aliquot.

- 8 -

CLAIMS

What is claimed is:

- 1 1. A method for detecting the presence of contamination in a nucleic acid amplification
2 reaction conducted on a sample, comprising the steps of:
3 conducting a first nucleic acid amplification reaction in said sample, wherein at least one
4 first nucleic acid primer used in said first nucleic acid amplification reaction comprises a first
5 portion that is complementary to a nucleic acid sequence in said sample, the amplification of
6 which is desired, and a second portion that is not complementary to said nucleic acid sequence;
7 conducting a second nucleic acid amplification reaction in said sample wherein at least one
8 second primer used in said second nucleic acid amplification reaction is complementary to said
9 second portion; and
10 detecting contamination in said sample as the presence of amplicon in said second nucleic
11 acid amplification reaction.
- 1 2. The method of claim 1, wherein said second portion is not complementary to any
2 contiguous nucleic acid present in said sample prior to said first nucleic acid amplification
3 reaction.
- 1 3. The method of claim 1, wherein said first nucleic acid amplification reaction is selected
2 from the group consisting of PCR, Q-PCR, and reverse-transcriptase PCR.
- 1 4. The method of claim 3, wherein said second nucleic acid amplification reaction is selected
2 from the group consisting of PCR, Q-PCR, and reverse-transcriptase PCR.
- 1 5. The method of claim 1, wherein said amplicon is detected by sequence-specific nucleic
2 acid probe capture.
- 1 6. The method of claim 1, wherein said first and second nucleic acid amplification reactions
2 are conducted simultaneously.

- 9 -

1 7. The method of claim 1, wherein said first and second nucleic acid amplification reactions
2 are conducted on DNA isolated from said sample.

1 8. A method for detecting contamination in a nucleic acid amplification reaction conducted
2 on a sample, comprising the steps of:

3 conducting a first nucleic acid amplification reaction in said sample using at least one
4 chimeric primer comprising a template-specific sequence and a 5' contamination detection
5 sequence;

6 conducting a second nucleic acid amplification reaction in said sample using at least one
7 primer that is substantially complementary to said contamination detection sequence; and

8 detecting an amplicon produced in said second nucleic acid amplification reaction, the
9 presence of which being indicative of contamination in said sample.

1 9. The method of claim 8, wherein said first nucleic acid amplification reaction comprises
2 two chimeric primers.

1 10. The method of claim 8, wherein said second nucleic acid amplification reaction comprises
2 two primers that are complementary to said contamination detection sequence.

1 11. The method of claim 8, wherein said sample is a biological sample.

1 12. The method of claim 11, wherein said sample is a stool sample.

1 13. The method of claim 8, wherein said amplification reaction is a polymerase chain reaction.

1 14. The method of claim 13, wherein said polymerase chain reaction is Q-PCR.

1 15. The method of claim 13, wherein said polymerase chain reaction is reverse-transcriptase
2 PCR.

1/5

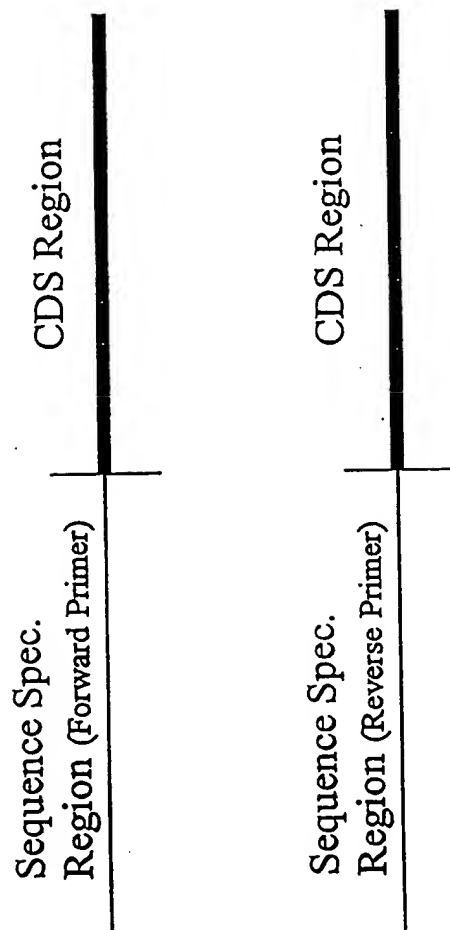


Fig. 1

2/5

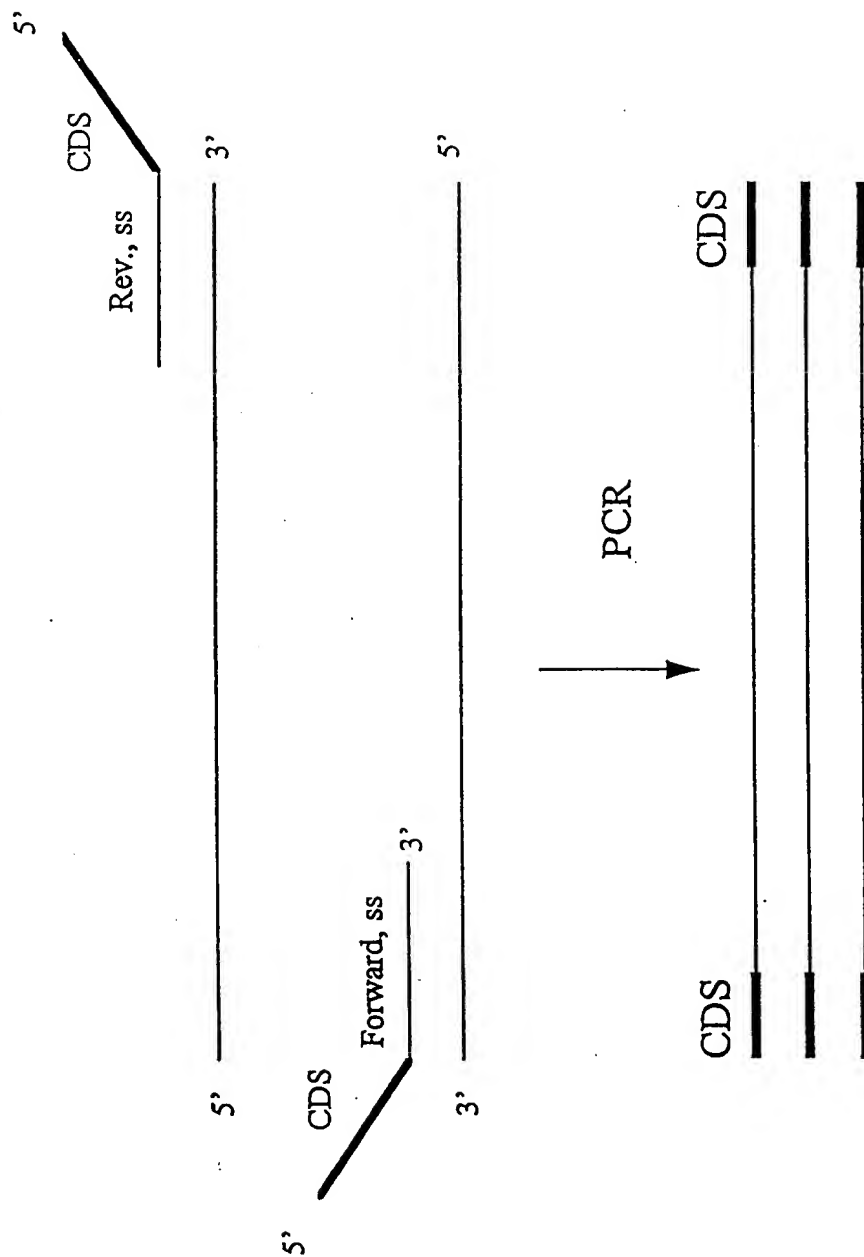


Fig. 2A

3/5

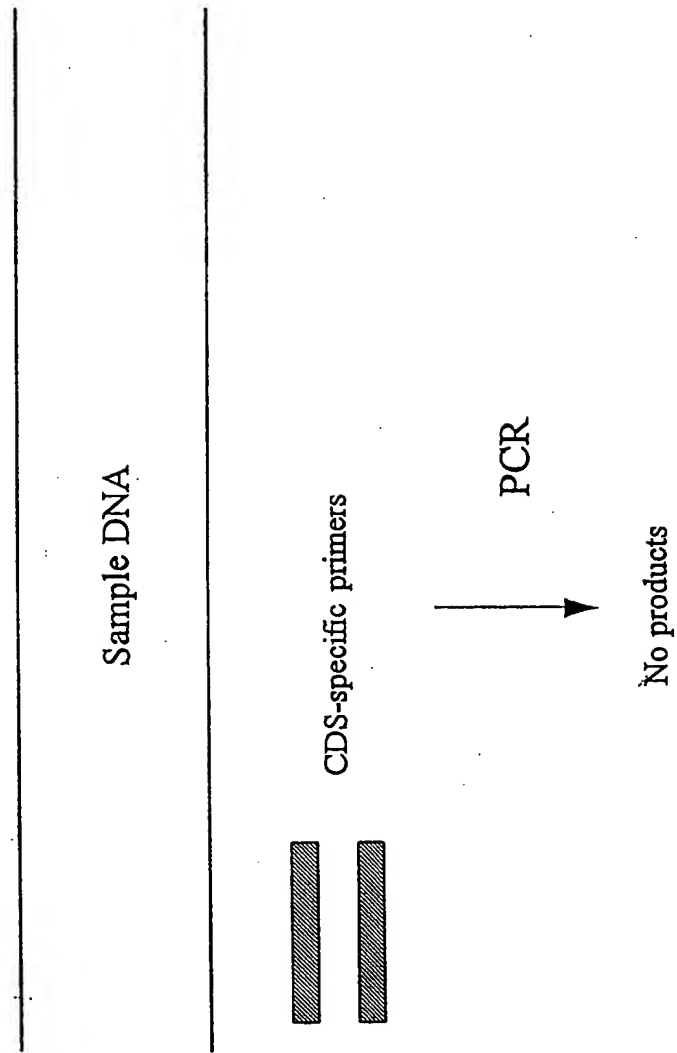


Fig. 2B

4/5

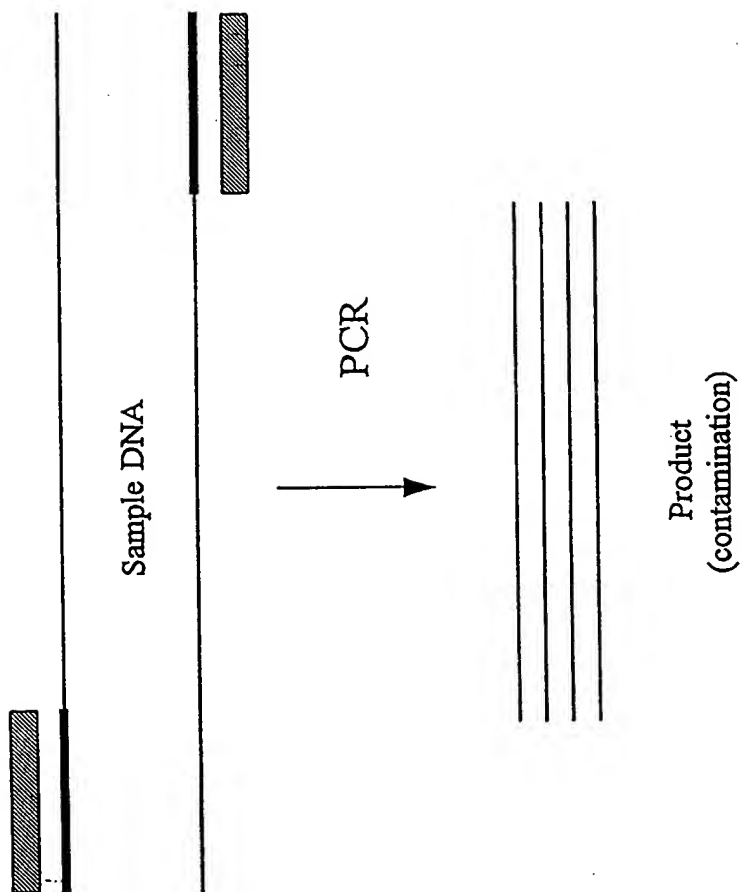
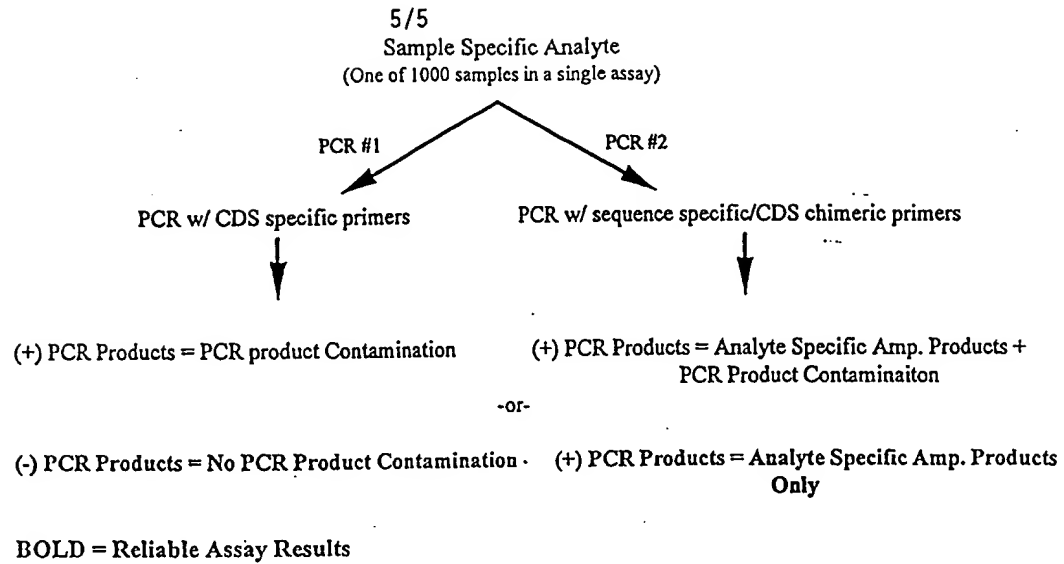


Fig. 2C

**FIGURE 3**

SEQUENCE LISTING

<110> Exact Laboratories, Inc.

<120> METHODS FOR DETECTING CONTAMINATION IN MOLECULAR
DIAGNOSTICS USING PCR

<130> EXT-010PC

<140>

<141>

<150> USSN 60/063,219

<151> 1997-10-23

<160> 3

<170> PatentIn Ver. 2.0

<210> 1

<211> 20

<212> DNA

<213> Artificial Sequence

<220>

<223> Description of Artificial Sequence:
Oligonucleotide Primer

<400> 1

gcggtcccaa aagggtcagt

20

<210> 2

<211> 28

<212> DNA

<213> Artificial Sequence

<220>

<223> Description of Artificial Sequence:
Oligonucleotide Primer

<400> 2

gattcctaca ggaagcaagt agtaattg

28

<210> 3

<211> 28

<212> DNA

<213> Artificial Sequence

<220>

<223> Description of Artificial Sequence:
Oligonucleotide Primer

<400> 3

taatggtgaa tatcttcaaa tgatttag

28

INTERNATIONAL SEARCH REPORT

International Application No
PCT/US 98/22406

A. CLASSIFICATION OF SUBJECT MATTER

IPC 6 C12Q1/68

According to International Patent Classification (IPC) or to both national classification and IPC

B. FIELDS SEARCHED

Minimum documentation searched (classification system followed by classification symbols)

IPC 6 C12Q

Documentation searched other than minimum documentation to the extent that such documents are included in the fields searched

Electronic data base consulted during the international search (name of data base and, where practical, search terms used)

C. DOCUMENTS CONSIDERED TO BE RELEVANT

Category *	Citation of document, with indication, where appropriate, of the relevant passages	Relevant to claim No.
X	WO 91 15601 A (US ARMY) 17 October 1991 see page 24, line 10-30; claim 1; figures 5,7 see page 13, line 10-14	1-15
A	WEIGHARDT F ET AL: "A SIMPLE PROCEDURE FOR ENHANCING PCR SPECIFICITY" PCR METHODS AND APPLICATIONS, vol. 3, no. 1, August 1993, pages 77-81, XP000609125 see the whole document	1-15
A	WO 91 17270 A (AMGEN INC) 14 November 1991 see abstract; figures 2,15 see page 7, line 1-18 see page 16, line 22 - page 23, line 10 see page 26, line 14 - page 27, line 18	1-15
-/--		

☒ Further documents are listed in the continuation of box C.

☒ Patent family members are listed in annex.

* Special categories of cited documents :

- "A" document defining the general state of the art which is not considered to be of particular relevance
- "E" earlier document but published on or after the international filing date
- "L" document which may throw doubts on priority claim(s) or which is cited to establish the publication date of another citation or other special reason (as specified)
- "O" document referring to an oral disclosure, use, exhibition or other means
- "P" document published prior to the international filing date but later than the priority date claimed

- "T" later document published after the international filing date or priority date and not in conflict with the application but cited to understand the principle or theory underlying the invention
- "X" document of particular relevance; the claimed invention cannot be considered novel or cannot be considered to involve an inventive step when the document is taken alone
- "Y" document of particular relevance; the claimed invention cannot be considered to involve an inventive step when the document is combined with one or more other such documents, such combination being obvious to a person skilled in the art.
- "Δ" document member of the same patent family

Date of the actual completion of the international search

5 March 1999

Date of making of the international search report

15/03/1999

Name and mailing address of the ISA

European Patent Office, P.B. 5818 Patentlaan 2
NL - 2280 HV Rijswijk
Tel. (+31-70) 340-2040, Tx. 31 651 epo nl,
Fax: (+31-70) 340-3016

Authorized officer

Reuter, U

INTERNATIONAL SEARCH REPORT

International Application No
PCT/US 98/22406

C.(Continuation) DOCUMENTS CONSIDERED TO BE RELEVANT		
Category	Citation of document, with indication, where appropriate, of the relevant passages	Relevant to claim No.
A	<p>EP 0 415 755 A (LIFE TECHNOLOGIES INC) 6 March 1991 see abstract see column 6, line 54 - column 9, line 5 -----</p>	1-15

INTERNATIONAL SEARCH REPORT

information on patent family members

Int. l. Application No

PCT/US 98/22406

Patent document cited in search report	Publication date	Patent family member(s)	Publication date
WO 9115601 A	17-10-1991	NONE	
WO 9117270 A	14-11-1991	AU 636499 B	29-04-1993
		AU 7987191 A	27-11-1991
		CA 2063432 A	02-11-1991
		CN 1059562 A	18-03-1992
		EP 0481069 A	22-04-1992
		FI 98932 B	30-05-1997
		IL 97981 A	18-06-1996
		JP 4506906 T	03-12-1992
		PT 97534 A,B	30-09-1993
		US 5427929 A	27-06-1995
		US 5650302 A	22-07-1997
EP 0415755 A	06-03-1991	AT 131878 T	15-01-1996
		CY 2073 A	11-09-1998
		DE 69024286 D	01-02-1996
		DE 69024286 T	30-05-1996
		DK 415755 T	22-01-1996
		ES 2080807 T	16-02-1996
		GR 3019092 T	31-05-1996
		HK 1000380 A	13-03-1998
		JP 2059842 C	10-06-1996
		JP 3091484 A	17-04-1991
		JP 7089932 B	04-10-1995
		US 5683896 A	04-11-1997